


MYOBRACE: MEASUREMENT SIZE CHART AND RULER GUIDE

Myobrace Size Conversion Chart: For use in place of standard Myobrace Ruler

MYOBRACE® SIZE AND CONVERSION CHART									
	Myobrace® Size	2-2 Upper (mm)	2-2 Lower (mm)	Canine Width Upper (mm)	Canine Width Lower (mm)	Arch Width 6-6 (mm)	Arch Depth A-P (mm)	Occluso-Guide™ Size™	Multi-P™ Size
SMALL	1	30.0	22.0	7.0	6.5	55.0	42.0	2½ - 3	30
	2	31.5	23.0	7.3	6.8	55.0	42.0	3 - 3½	35
MEDIUM	3	32.5	23.8	7.5	7.0	60.0	48.0	4 - 4½	40
	4	33.5	24.5	7.7	7.2	60.0	48.0	4½ - 5	45
	5	35.0	25.5	8.0	7.5	60.0	48.0	5 - 5½	50
LARGE	6	36.5	26.7	8.5	8.0	60.0	48.0	6 - 6½	55
	7	38.0	28.0	9.2	9.0	66.0	52.0	7 - 7½	60

Measurement Guidelines

When measuring, you will need to measure from distal-of-lateral to distal-of-lateral (mid-labial) #7-10, or the four upper anterior teeth. If the patient is missing their laterals you will want to measure from mesial of k-9 to mesial of k-9. Include all crowding and spacing. If the measurement is: <32mm it will be a SMALL, 33-36mm a MEDIUM, and >37mm will be a LARGE. If the patient measures between sizes select the larger size.

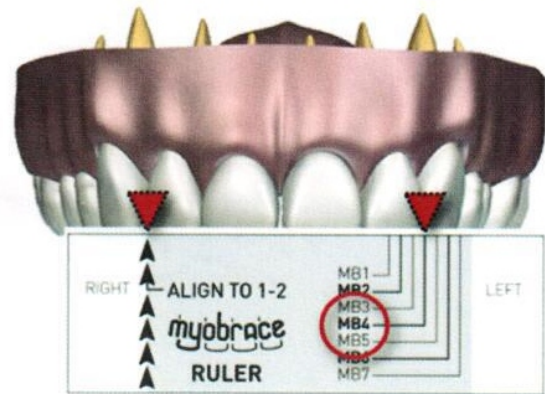

Myobrace Rulers for Office Use and Lamination

Use with Myobrace T Series

